

CH05 Chiller and AFA ultraCUBE Water Hose Replacement Warning



If this hose kit is part of an upgrade of existing CH03 or CH04 water chiller or WCS system the black rubber hoses previously used must be discarded and replaced with beige hoses from this kit.



CH03 Chiller




CH04 Chiller




CH05 Chiller

Existing hose is not adequately rated for the water pressure generated by the CH05 or AFA ultraCUBE.

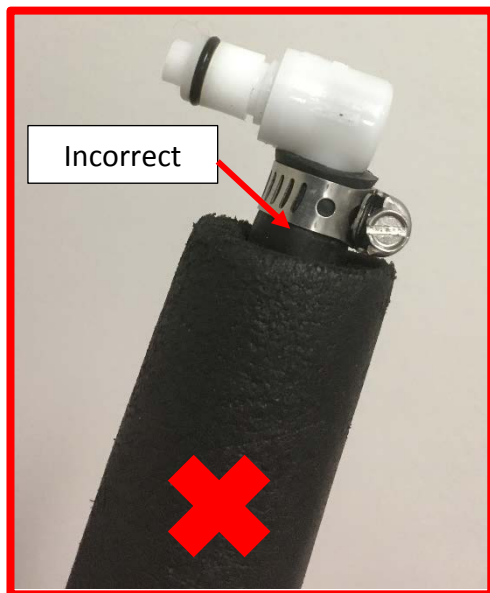

Remove and discard these hoses.

**Includes part numbers:
190321, 500102, 500173, 500174**

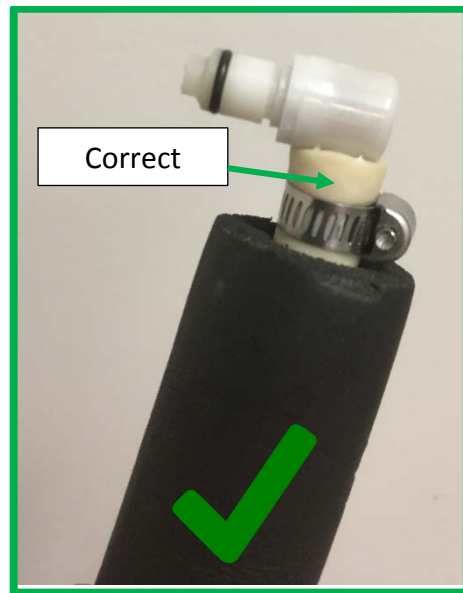

Install new hoses from this kit.

**Includes part numbers:
500575, 500576**